

Upholstery Maintenance Guide



Cleaning

A-dec recommends cleaning the upholstery using a solution of mild, non-ionic detergent and water, or commercially available cleaners containing no alcohol, bleach, or ammonia. Dishwashing detergent is usually non-ionic. Mix just enough detergent to allow for good cleaning without leaving a soapy film on the surface. Never use abrasive cleansers, scrubbing pads, or other abrasive applicators. They can permanently scratch or otherwise damage upholstery surfaces.

Barrier Protection

There is low risk that pathogens will be transferred from chair upholstery. However, since interpretation of regulatory agency guidelines commonly implies some form of infection control related to upholstery, A-dec recommends barrier protection instead of relying on chemicals. Barriers significantly extend the life of your chair upholstery and will help to preserve its luxurious look and soft feel. In the USA, barriers must be produced under the Current Good Manufacturing Practice (CGMP) as specified by the U.S. Food and Drug Administration (US FDA). For regions outside the USA, refer to the medical device regulations specific to your location.



CAUTION If you use barriers, always replace the barrier film after each patient.

Disinfectants



CAUTION Minimize the use of surface disinfectants on the upholstery. Chemicals can cause deterioration and shorten the life of upholstered surfaces.

If a barrier becomes compromised or if visible contamination is present on an upholstered surface, an EPA-registered hospital grade disinfectant should be used. Disinfectants may cause damage and no disinfectant should be considered harmless to dental equipment. The degree of damage is influenced by several factors including, but not limited to: the chemical(s) used, level of exposure, and length of time. Your dental dealer has knowledge and experience with barriers and disinfectant products and will provide you with good advice.

Upholstery Life and External Factors

A-dec uses industry standard and other proprietary tests to ensure robust upholstery design. We also use this testing to validate how the upholstery will react and wear over time when exposed to a variety of external factors.

These factors can impact the longevity of upholstery and cause premature deterioration of its surfaces. It is the responsibility of the clinician to understand and regulate exposure to these elements to ensure the maximum useful life of A-dec upholstery. The most common factors include:


- chemicals (disinfecting/cleaning)
- ultraviolet light (sunlight)
- humidity, heat, and high mineral content water
- equipment lubricants, materials used in dental procedures, and sharp objects
- abrasive applicators and cleaners
- body oils and non-color-safe materials

Deterioration


The common indicators of the external factors affecting upholstery life include: changes in the texture, hardening, cracking, peeling, splitting, and discoloration of the upholstered surface.




Cracking

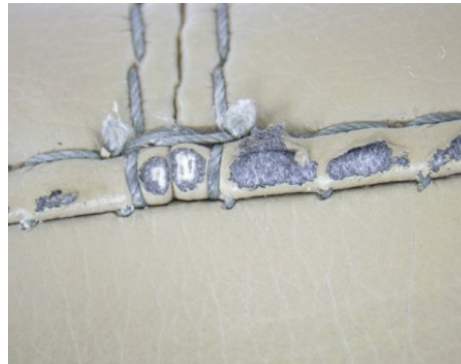

Peeling

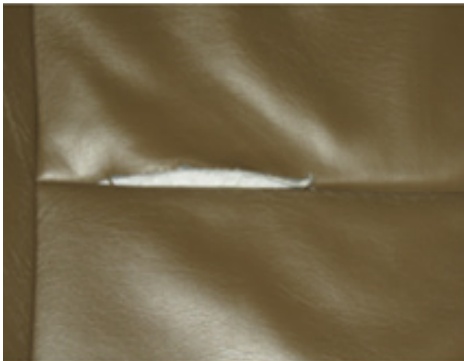

Tearing at the Seam

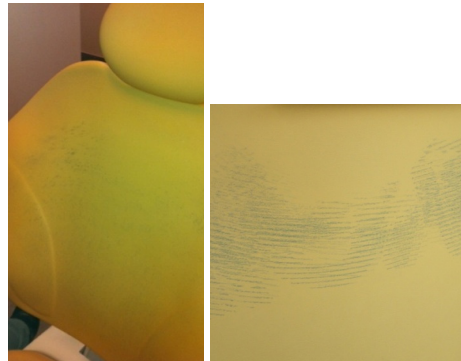

Chemical Staining / Dye Transfer

Texture Deterioration

Chemical disinfection is the largest contributor to upholstery deterioration. Although symptoms can vary, seamless upholstery may become glossy in appearance, the texture may become smooth, and the upholstery may harden and crack over time. Sewn upholstery may become tacky and dull, eventually cracking, tearing at the seam, or peeling, which exposes the backing material.

Extreme heat and humidity can also be detrimental to the construction of A-dec upholstery over time. Although the upholstery is engineered to withstand most environmental conditions, prolonged extreme exposure can break down the upholstery substrate.

Abrasions and Punctures

Abrasive applicators or abrasive cleaning products can damage the upholstery by scratching and roughing up the textured surface. Sharp objects can puncture and tear the material.

Discoloration

Exposure to sunlight over extended periods of time can fade the upholstery. The use of curtains or blinds, when appropriate for window facing chairs, will help limit exposure to ultraviolet light (UV).

Chemical Staining and Dye Transfer

Some chemicals will discolor and stain the upholstery. Non-color-safe materials, combined with heat and perspiration, can transfer dye onto the upholstery.

Chemical Disinfectant Testing on A-dec Upholstery

A-dec continually tests (based on ASTM and ISO standard test methods) for interactions between chemical disinfectants and A-dec equipment surfaces. These standard test methods specify test conditions to reveal the effects that may occur with repeated applications of a surface disinfectant over time. A-dec conducts this testing to help select the best suited materials for chemical resistance, as well as other important characteristics. In the clinical environment, the manual process of surface cleaning and disinfection is subject to significant variation, which may exceed the manufacturer's recommended instructions. At A-dec, our upholstery testing methods incorporate a reasonable amount of this exposure, as well as other factors (e.g., abrasion, UV, etc.).

We recommend using the active agent as an indicator for chemical aggressiveness.

In our experience, surface disinfectants containing any of the following active agents are among the **most** detrimental to the life of the upholstery:

- isopropanol or isopropyl alcohol (greater than 25% by volume)
- hydrogen peroxide
- chlorine compounds (such as sodium hypochlorite)

Our findings indicate that the following active agents are among the **least** detrimental (but may not be harmless) to the life of the upholstery:

- high-dilution water-based phenolic compounds
- quaternary ammonium compounds ("quats" with less than 25% alcohol by volume)

For more information, please contact A-dec Customer Service or see the *A-dec Equipment Asepsis Guide* (p/n 85.0696.00), which is available at www.a-dec.com and ships with all A-dec products.



A-dec Headquarters

2601 Crestview Drive
Newberg, Oregon 97132
United States
Tel: 1.800.547.1883 within USA/CAN
Tel: +1.503.538.7478 outside USA/CAN
www.a-dec.com

International Distribution Centers

A-dec United Kingdom

Austin House
11 Liberty Way
Nuneaton, Warwickshire CV11 6RZ
England
Tel: 0800.ADEC.UK (2332.85) within UK
Tel: +44.(0).24.7635.0901 outside UK

A-dec Australia

Unit 8
5-9 Ricketty Street
Mascot, NSW 2020
Australia
Tel: 1.800.225.010 within AUS
Tel: +61.(0).2.8332.4000 outside AUS

86.0501.00 Rev C
Date of Issue 2020-02-12
Copyright 2020 A-dec Inc.
All rights reserved.