


A-dec's Simple Solution To Drying Preps.

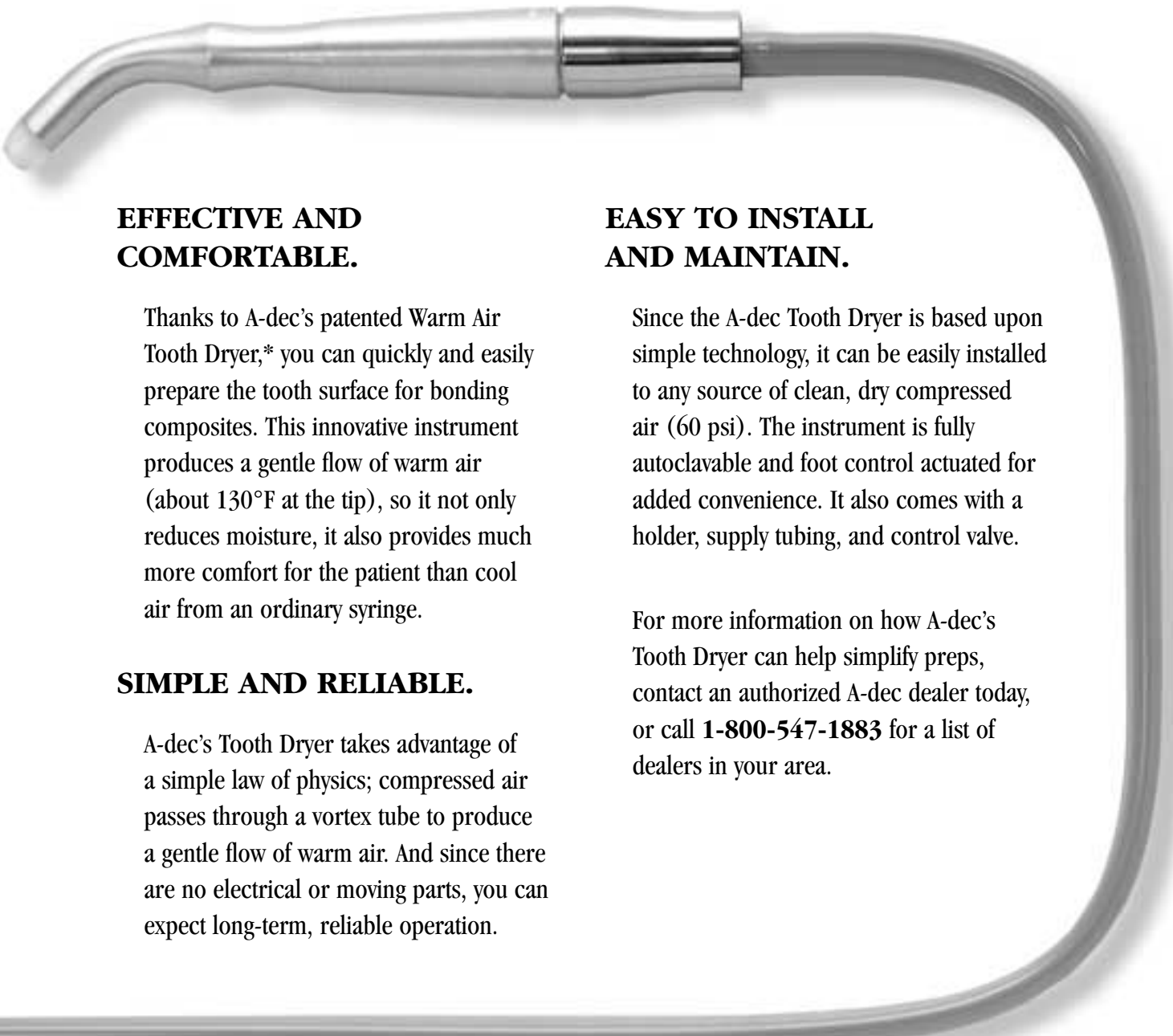

EFFECTIVE AND COMFORTABLE.

Thanks to A-dec's patented Warm Air Tooth Dryer,* you can quickly and easily prepare the tooth surface for bonding composites. This innovative instrument produces a gentle flow of warm air (about 130°F at the tip), so it not only reduces moisture, it also provides much more comfort for the patient than cool air from an ordinary syringe.

SIMPLE AND RELIABLE.

A-dec's Tooth Dryer takes advantage of a simple law of physics; compressed air passes through a vortex tube to produce a gentle flow of warm air. And since there are no electrical or moving parts, you can expect long-term, reliable operation.

EASY TO INSTALL AND MAINTAIN.

Since the A-dec Tooth Dryer is based upon simple technology, it can be easily installed to any source of clean, dry compressed air (60 psi). The instrument is fully autoclavable and foot control actuated for added convenience. It also comes with a holder, supply tubing, and control valve.

For more information on how A-dec's Tooth Dryer can help simplify preps, contact an authorized A-dec dealer today, or call **1-800-547-1883** for a list of dealers in your area.



A-dec Inc.
PO Box 111, 2601 Crestview Drive
Newberg, OR 97132 USA
Phone: 1-800-547-1883, Fax: 503-538-0276

Outside USA/Canada, contact:
A-dec International Inc.
Phone: 503-538-7478, Fax: 503-538-5911
©A-dec Inc., 1999. All rights reserved. 85.6001.00/1mx/2.5M/3-99

*US Patent No. 4,973,246

